MANUAL DE INSTRUCCIONES / USER GUIDE

CAMILLA CUCHARA / SCOOP STRETCHER

DT05_01_07 // REV. 2026/05

00 ÍNDICE / INDEX

REV. 2026/05

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01.2. Using this guide

This manual provides the instructions for use and maintenance of the product, together with technical information, operating procedures, spare parts and safety information.

Before placing the product into service, it is recommended that this manual be read carefully, regardless of the level of experience previously acquired with similar devices, in order to prevent damage caused by improper use.

Do not lose this document. It must remain available for any query that may arise from clinical personnel. If the instructions for use refer to a product different from the one received, the manufacturer must be contacted immediately before using the device.

Remember: proper use and maintenance are essential for the correct operation of the device.

01.3. Legend of symbols

SYMBOL	EXPLANATION / DESCRIPTION
	MANUFACTURER symbol. This symbol is accompanied by the name and address of the manufacturer, adjacent to the symbol (PRODUCTOS METÁLICOS DEL BAGÉS S.L., Ctra. C-16 Km 59.5, 08650 Sallent (Barcelona)).
	Indicates the manufacturer's reference number to identify the medical device. PROMEBA, S.L. uses this symbol to set each internal reference for each configuration and business variant
	Indicates the manufacturer's serial number to identify a specific medical device.
	Indicates the manufacturing date. The symbol must be accompanied by a manufacturing date (yyyy-mm), adjacent to the symbol.
	It is placed to inform that the product is a "Medical Device"
	CE symbol without the intervention of a Notified Organism, as a Medical Device classified as Class I according to the EU Regulation 2017/745 on Medical Devices.
	Symbol for the Unique Device Identifier.
	Symbol "See instructions for use or operating instructions"
	Symbol "caution". This symbol is placed to warn of the need for the user to refer to important precautionary information in the operating instructions, such as warnings and cautions not otherwise found on the label.

01 INTRODUCTION


REV. 2026/05

Symbol "Caution". For a general warning.

Do not use in MRI's

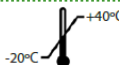Indicates the temperature limits. The upper and lower temperature limits should be indicated adjacent to the horizontal lines.

01.4. Customer service

For information regarding the correct interpretation of the instruction manual, use, maintenance, installation and restoration of the product, please contact Promebea customer service at T. (+34) 93 837 12 00, mobile (+34) 608 15 31 93, e-mail promebea@promebea.com, or write to PROMEBEA, S.L. — Ctra. C-16 Km 59.5 · 08650 Sallent (Barcelona) · SPAIN.

Whenever any matter relating to a product is reported, the product reference (**REF**) and the serial number (**SN**) will be requested. It is strongly recommended to keep the label containing this information, which is provided in the original packaging of the product.

01.5. Service life

If used as described in the following instructions, this device has a useful life of 8 years from the date of purchase.

Promebea, S.L. accepts no liability for malfunction or damage caused by the use of devices that have exceeded the maximum permitted service life.

01.6. Disposal

When the devices are no longer fit for use, and provided they have not been contaminated by any specific agent, they may be disposed of as normal solid waste. Otherwise, the disposal regulations in force must be followed.

01.7. Reference standards

REFERENCE	TITLE OF DOCUMENT
EU Regulation 2017/745	EU Regulation on Medical Devices
UNE-EN 1865-1:2012	Patient handling equipment used in road ambulances - Part 1: General stretcher systems and patient handling equipment


As a distributor or end user of the products manufactured and/or distributed by Promebea, S.L., you are strictly obliged to know the legal provisions in force in the country of destination of the goods, applicable to the devices to be supplied (including the regulations relating to technical specifications and/or safety requirements) and, therefore, understand the necessary requirements to ensure compliance of the products themselves with all the legal requirements of the territory.

01 INTRODUCTION

REV. 2026/05

01.8. Labelling

Each product is fitted with an identification label, placed on the device itself and/or on the original packaging. It must never be removed or covered. This label includes the data of the company which produced the device () the product code / reference () of the product, the manufacturing date () the serial number () the UDI-PI code () and the CE distinctive (). Keep this label so that you can inform the distributor when required.

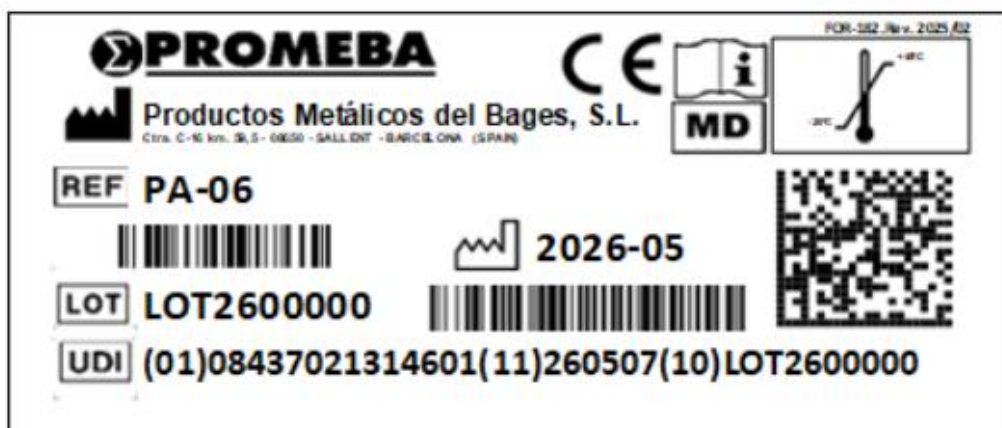The UDI-PI code is made up of the following fields:

- (01)84337021314601: product UDI code.
- (11)260507: information about the batch release date in the format *yyymmdd*, where *yy* is the year of release, *mm* is the month and *dd* is the day.
- (10)LOT2600000: (*LOT*) identifies the batch number, providing information about an internal document assigned by the IT system: year (*yy*) and consecutive number (*nnnnn*).

The original packaging of the product is labelled as shown below. It contains the same information as the product label, together with the product description in Spanish, English, Italian, French and German.

01 INTRODUCTION

REV. 2026/05

Σ PROMEBA

REF PA-06 

 2026-05

LOT LOT2600000 

MD   

TABLERO ESPINAL PA-06 X-RAY
X-RAY SPINE BOARD PA-06
TABLE DE SPINE PA-06 X-RAY
RÜCKNBRETT PA-06 X-RAY
TABLA DE SPINA PA-06 X-RAY

UDI (01)08437021314601(11)260507(10)LOT2600000

 Productos Metálicos del Bages, S.L.
Ctra. C-16 km. 59,5
08650 - SALLENT
BARCELONA
(SPAIN)

FCR-18.3.Rev. 2025/02

02 INTENDED USE

REV. 2026/05

02.1. Intended use

The SCOOP STRETCHER family of products is designed to lift and load patients onto immobilization and transport devices such as ambulance cots. With the use of the device, the patient can be picked up without changing position, and the trunk and limbs can be kept aligned during loading.

02.2. Contraindications and side effects

Absolute or strongly restricted contraindications:

The scoop stretcher is contraindicated in patients with isolated penetrating trauma (gunshot or stab wounds) without neurological deficit, since its application delays transfer to the hospital without providing any proven clinical benefit and may worsen the prognosis (NAEMSP/ACS-COT).

It must not be used as a definitive transport device or during prolonged transfers, especially on uneven terrain, since its rigidity transmits vibrations and impacts to the patient. For transport, the patient must be transferred to a main stretcher or, preferably, to a vacuum mattress.

It must not be used for vertical lifting, rope rescue or high-angle rescue: the scoop stretcher is not designed as a lifting device with slings, cranes or hoisting systems, unless expressly indicated by the manufacturer.

It must not be used with patients whose weight exceeds the maximum load capacity stated by the manufacturer in the technical specifications, or whose anthropometry exceeds the dimensions of the device.

Relative contraindications:

Patients with severe kyphosis, ankylosing spondylitis, advanced scoliosis or marked thoracic deformities, in whom forcing the supine position on a rigid surface may cause iatrogenic injuries, including spinal cord injury. In these cases the patient must be immobilized in the position of maximum comfort, with adequate padding.

Pregnancy from the 20th week of gestation onwards, due to the risk of aortocaval compression in the supine position. If its use is essential, the stretcher must be tilted 15-30° to the left side (left uterine displacement) or the right side must be wedged.

Morbid obesity, elderly patients with severe osteoporosis, cachexia or skin fragility, due to the increased risk of tissue pinching when closing the two halves of the device and of early development of pressure injuries.

Combative or agitated patients, or patients with altered level of consciousness (intoxication by alcohol or drugs, traumatic brain injury), in whom the struggle against the straps may cause more harm than the potential benefit.

Severe respiratory failure or risk of bronchoaspiration, since the forced supine position may further compromise ventilation and hinder airway protection. In these patients, lateralization or elevation of the head end must be considered.

02 INTENDED USE

REV. 2026/05

Patients with suspected unstable pelvic fracture when the application manoeuvre might imply inappropriate lateral mobilization; in such cases, the technique must be performed with minimal movement ("skin to scoop").

Pediatric patients whose height is less than the minimum adjustable length of the stretcher, in whom full support of the head, trunk and lower extremities cannot be guaranteed.

Adverse effects:

The adverse effects identified from the clinical evaluation of the SCOOP STRETCHER product family are the following:

- Pain or discomfort caused by prolonged contact with the rigid surface.
- Impaired respiratory function caused by immobilization or by excessive strap tension.
- Pressure sores (decubitus ulcers) on bony prominences when application time exceeds the strictly necessary period.
- Pinching or entrapment of the patient's skin and soft tissues between the two halves of the device during closure, particularly on lateral aspects (flanks, buttocks).
- Hypothermia or thermal discomfort due to contact with the metal (aluminium) structure at low ambient temperatures.
- Need for additional radiological examinations to rule out imaging artefacts caused by the metal structure of the device.
- Increased pain in patients with lumbar or pelvic injuries due to the lack of continuous support in the central area, characteristic of the "scoop" design.
- Risk of patient drop in the event of incorrect engagement of the upper and lower latches.

02.3. Physical requirements of the operators

The Promeba SCOOP STRETCHER product family is intended exclusively for professional use. The device is expected to be used by rescue personnel with knowledge of the handling and transport of polytrauma patients.

Operators must be trained to transport patients in an efficient, effective and safe manner, and must meet the following minimum requirements:

- Physical capacity for operating the device.
- Strong back, arms and legs to lift and carry the device.
- Good muscular coordination.
- Capacity to grip the device firmly with both hands.

It is recommended the employment of operators with strength, balance, coordination and common sense.

Patient loading procedures for extremely heavy patients, operations in rough terrain and in particular situations may require more operators.

The capacities of the various operators must be considered before determining his role in the employment of the board.

Untrained users may injure themselves and cause harm to themselves or the patient. Use of the device should be entrusted exclusively to trained and qualified personnel.

03 WARNINGS

REV. 2026/05

03.1. General warnings

1. The product must be used only by trained personnel who have received specific training on this device and not on similar products.
2. Establish a maintenance and periodic-inspection programme and identify a person in charge. The person responsible for the routine maintenance of the device must comply with the basic requirements provided by the manufacturer in these instructions for use. It is recommended that this person check at least every 6 months whether updated user manuals are available or whether any changes have been made to the product. This information is freely available on the website <http://promeba.com/>.
3. All maintenance activities must be recorded and documented with the corresponding technical intervention reports. The documentation must be retained for at least 10 years from the end of the device's service life and must be made available to the competent authorities and/or the manufacturer upon request.
4. Promeba, S.L. is always at your disposal to plan product training sessions.
5. Before performing any operation on the product (training, installation, use), the operator must carefully read the enclosed instructions, paying particular attention to the safety precautions and to the procedures to be followed for correct installation and use.
6. If the instructions belong to a device other than the one received, notify the manufacturer immediately and refrain from using the device.
7. In case of any doubts regarding the correct interpretation of the instructions, contact Promeba, S.L. for the necessary clarification.
8. Do not allow untrained persons to assist during the use of the device, as they could cause injury to the patient or to themselves.
9. Before each use, the perfect operating condition of the device must be verified as specified in the instruction manual. If any damage or anomaly that could in any way influence the correct functioning and the safety of the device, the patient and/or the user is detected, the device must be immediately removed from service and the manufacturer must be contacted.
10. If any failure or malfunction of the device is detected, it must be immediately replaced with a similar item so that immobilization and rescue procedures are guaranteed without interruption.
11. Use of the device in any way other than that described in this manual is prohibited.
12. Do not alter or modify the product in any way; any such interference may cause malfunction and injury to the patient and/or the operator.
13. The product must not be tampered with in any way (modification, adjustment, addition, replacement). In such cases, all liability for any malfunction or injury caused by the product itself shall be disclaimed; in addition, the CE marking and the product warranty shall be deemed null and void.
14. Ensure that all necessary precautions are taken to avoid hazards arising from contact with blood or bodily fluids.

03 WARNINGS

REV. 2026/05

15. Record and store with these instructions: batch number, place and date of purchase, first date of use, date of inspections, names of users and any comments.
16. During use of the device, the assistance of qualified personnel must be guaranteed.
17. Position and adjust the device with care, so as not to obstruct the rescuers or any other rescue equipment.
18. The product must not be exposed to, or come into contact with, thermal sources of combustion or flammable agents; it must be stored in a dry, cool place protected from light and sunlight.
19. Store and transport the device in its original packaging.
20. Do not store the product beneath other materials that may damage the product structure due to their weight.
21. Attention: laboratory tests, post-production tests and instruction manuals cannot always foresee every possible use scenario. This means that in some cases the performance of the product may differ significantly from the results obtained to date. The instructions are continuously updated and are subject to strict supervision by fully qualified personnel with the appropriate technical training.
22. As a distributor or end user of products manufactured and/or marketed by Promeba, S.L., it is strictly required that you have basic knowledge of any legal requirements in force in the country of final destination of the goods applicable to the devices included in this supply (including laws and standards relating to technical specifications and/or safety requirements). You are therefore also strictly required to have the knowledge necessary to guarantee all aspects relating to the full conformity of the products with the regulations of the relevant territory.
23. Promptly notify Promeba, S.L. regarding any revisions that must be carried out by the manufacturer in order to guarantee the conformity of the product with the legal specifications of the territory (including those resulting from rules and/or regulatory provisions of another nature).
24. Act with due care and diligence, and contribute to ensuring the conformity with the general safety requirements of all devices placed on the market in the territory, providing end users with all the necessary information for the periodic inspection of their devices, as specified in the user manual.
25. Actively contribute to the safety controls of products sold by communicating any relevant risk-analysis information to both the manufacturer and the competent authorities so that the necessary actions can be taken promptly.
26. Be aware that, in the event of failure to comply with the above requirements, you shall be deemed fully responsible for all damages that may occur. We therefore expressly disclaim any liability and/or obligation for non-compliance with the present 'Regulatory provisions'.

03 WARNINGS

REV. 2026/05

03.2. Specific warnings

1. Use only accessories/spare parts that are original or approved by Promeba, S.L. when performing any operation, so as not to cause any alteration or modification to the device. Otherwise, no liability shall be assumed for correct operation or for any resulting damage to the patient or the operator, and the warranty shall be deemed null and void, as well as the conformity with Regulation (EU) 2017/745.
2. Always respect the maximum capacity of the device, as indicated in this user manual. Maximum load capacity means the total weight distributed according to human anatomy. When determining the total weight load on the product, the operator must take into account the weight of the patient, equipment and accessories. In addition, the operator must consider that the overall dimensions of the patient should not reduce the functionality of the device.
3. The device and all its components, after washing, must be allowed to dry completely before storage. Do not use tumble dryers to dry the device.
4. Avoid contact with sharp objects.
5. Do not use the device if it shows cracks, deformations or abrasions on the aluminium structure, or excessive wear of the latches, hinges, extension joints or straps.
6. Before lifting, make sure that operators have a firm grip on the device.
7. Avoid dragging the device over rough surfaces.
8. Do not use as a parking device; transfer the patient to a vacuum mattress or to the primary ambulance stretcher as soon as possible.
9. A minimum of two operators in adequate physical condition is required to use the device. They must have strength, balance, coordination and common sense, and must be trained in the correct operation of the device.
10. For the loading of particularly heavy patients and for rescue operations on steep terrain or in unusual circumstances, the presence of additional operators is recommended.
11. Do not operate if the weight has not been distributed properly.
12. Pay close attention to possible obstacles (water, ice, debris, etc.) along the route, as they could cause the operator to lose balance and compromise the proper functioning of the device. If the path cannot be cleared of obstacles, choose an alternative route.
13. Operators must work together and maintain control of the device throughout the entire operation.
14. Follow the procedures and standard protocols for handling the patient in an emergency.
15. Lift only the weight you can safely support.
16. Never leave the patient on the device unattended; injury could occur.

03 WARNINGS

REV. 2026/05

17. Use of the device without proper patient immobilization may cause serious harm. Always ensure that the patient is properly immobilized with the straps supplied before using the device.
18. The rescuer must assess the actual need to apply this type of device in accordance with local guidelines.
19. Perform rescue simulations with a stretcher and a load simulating a patient and accessories before placing the device into service.
20. The device is intended to be in contact with the patient through the patient's clothing. In the event of direct contact with the skin, place a protective surgical drape in between to prevent contamination.
21. During lifting, the stretcher must be kept horizontal with respect to the ground; any abnormal tilt may cause serious harm to the patient, the device and the operator.
22. To preserve the service life of the device, protect it from UV radiation and adverse weather conditions.
23. Do not machine-wash the device.
24. Do not use solvents or stain removers.
25. Do not use bleach to disinfect; use an alcohol-based disinfectant.
26. Condensation, water, ice and dust may compromise the proper operation of the device.
27. Before closing the two halves, visually and manually check that the patient's skin, clothing, straps and life-support equipment are not trapped between the halves. Risk of skin or tissue injury.
28. Before any lifting manoeuvre, check by visual inspection and by performing a pull test that both latches have engaged correctly. Accidental release of one of the latches may cause the patient to fall.
29. Adjust the length of the stretcher to the patient's height before applying it. Insufficient length may leave the head or feet unsupported; excessive length may prevent the two halves from engaging correctly.
30. Do not use with cranes, slings, hoisting systems, rope rescue systems or in high-angle or vertical rescues. The handles and latches are not designed as external mechanical load points.
31. Before folding the stretcher, make sure that the extension joints are fully exposed. Folding it without doing so may damage the device irreparably.
32. At low ambient temperatures, the aluminium structure may transmit cold to the patient. Place a thermal blanket or sheet between the stretcher and the patient.

03.3. Residual risk

The residual risks listed below have been identified solely in relation to the intended use of the device and cannot be fully eliminated, even when all the warnings in this manual are

03 WARNINGS

REV. 2026/05

observed. They can only be minimized through proper operator training, strict compliance with the application procedure and limitation of the time the patient spends on the device:

1. Risk of cross-infection due to inadequate cleaning and disinfection procedures.
2. Risk of patient drop during lifting and manual transport, even when the patient is properly secured and the latches are correctly engaged.
3. Risk of accidental release of one of the latches (head or foot end) due to impacts, vibrations or improper closure.
4. Risk of pinching of the patient's skin, clothing or soft tissues between the two halves during closure, inherent to the scoop design.
5. Risk of unintended movement of the injured area during application and transfer, even when correctly performed by trained personnel, particularly in cases of suspected unstable pelvic fracture or spinal injury.
6. Risk of airway compromise in unconscious patients due to the supine position and immobilization with straps.
7. Risk of pain, respiratory impairment and early onset of pressure injuries associated with contact with the rigid surface, only minimizable by limiting application time to the strictly necessary period.
8. Risk of hypothermia or thermal discomfort due to direct patient contact with the aluminium structure at low ambient temperatures.
9. Risk of insufficient number of operators in real rescue conditions relative to the recommended optimum, which may compromise the safety of the patient and of the operators themselves. The appropriate number must be assessed according to local regulations and the characteristics of the intervention.

04 PRODUCT DESCRIPTION *REV. 2026/05*

04.1. Main components

Feature	10XLP	PA-05	PA-179
Image			
Trade name	CAMILLA CUCHARA 10XLP	CAMILLA CUCHARA ALUMINIO	CAMILLA PALAS Y TABLERO ESPINAL
Intended use	Adult	Adult	Adultt
Length (max)	2.090 mm	2.140 mm	1.860 mm
Length (min)	1.670 mm	1.670 mm	1.860 mm
Width	420 mm	420 mm	430 mm
Thickness (height)	70 mm	70 mm	45 mm
Folded dimensions	1.200x420x70 mm	1650x420x70 mm	
Weight	7,5 kg	8 kg	8,5 kg
Maximum load	200 kg	170 kg	170 kg
Material	PE	Aluminio	PE
Color	Yellow	Aluminium	Green
Straps included	3	3	3

05.1. Transport and storage

Before transporting the product, make sure it is correctly packaged. During transport, always secure the load. If stacking is necessary, follow the scheme indicated in the packaging documentation.

Transport the load levelled and in compliance with the general rules for the transport of loads, ensuring there is no risk of impacts or falls during the journey.

To unpack the device, place the box on a flat, stable surface and carefully open the seal. Keep the original packaging to use for any further transport and for storage.

The device must be stored in a cool, dry place, away from direct sunlight, and must not come into contact with substances or chemical agents that may damage it or reduce its safety performance. Storage temperature: from –20 °C to +40 °C.

Damage to the product caused during transport and handling is not covered by the warranty. Repair or replacement of damaged parts is the customer's responsibility.

The product must not be manipulated in any way (modification, adjustment, addition, substitution). In such cases, all responsibility for any malfunction or injury caused by the product itself will be disclaimed; furthermore, CE certification and product warranty will be considered void.

05.2. Preparation

On receipt of the product:

- Remove the packaging and display the material so that all components are visible.
- Check that all the components/pieces are present.

Before each use, the device must be inspected to detect any operating anomalies and/or damage caused by transport and/or storage. In particular, check:

- General functionality of the device.
- Cleanliness of the device (remember that the failure of cleaning may cause the risk of cross infections).
- Absence of cracks, deformations, abrasions or breaks in the aluminium structure and in the straps.
- Integrity of all components.
- Correct operation of the length-adjustment buttons and of the upper and lower latches.

If the above conditions are met, the device may be considered ready for use; otherwise, it must be immediately removed from service and the manufacturer must be contacted.

05.3. Functioning

Guidelines for the use of the scoop stretcher

Before using the scoop stretcher, carefully read this manual and the manuals of all accessories intended to be used with the device (restraint straps, spider straps, cervical collar, head immobilizer, vacuum mattress).

05 OPERATION

REV. 2026/05

Strictly follow the guidelines of your emergency service and the local or regional protocols before placing the patient on the stretcher. A primary medical assessment of the patient must be carried out before each intervention.

Carefully assess that the weight to be moved does not exceed the operators' capacity. If it does, additional personnel will be required for transport.

Assess the actual need to use the scoop stretcher according to local guidelines, as well as the suitability of applying a cervical collar, head immobilizer or other means of immobilization in addition to the straps.


Length adjustment:


Open the stretcher and place it on the floor so that the two halves are perfectly aligned, with the upper and lower latches engaged (not separated).

1. Unlock the adjustment button on both sides of the stretcher and hold it in the vertical position.

2. Extend or shorten the length of the stretcher by sliding the profile to one of the available positions, according to the patient's height.

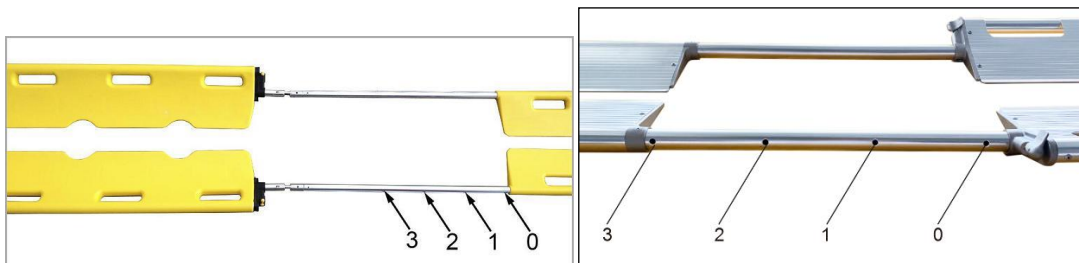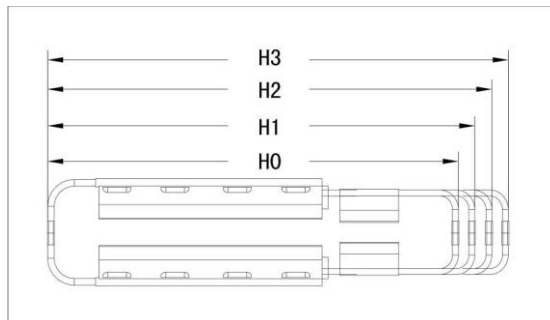3. Release the button on both sides so that it returns to its locked position. Check that the profile is firmly fixed in the selected position.


Insufficient length may leave the patient's head or feet unsupported; excessive length may prevent the two halves from engaging correctly when re-joined.

The maximum adjustment length of the stretcher is 200cm, that is, the length cannot be modified by more than 200cm.

Separation of the two halves:

Simultaneously press the buttons of the upper (head) and lower (foot) latches and separate the two halves of the stretcher. Place each half on either side of the patient, parallel to the longitudinal axis of the body.

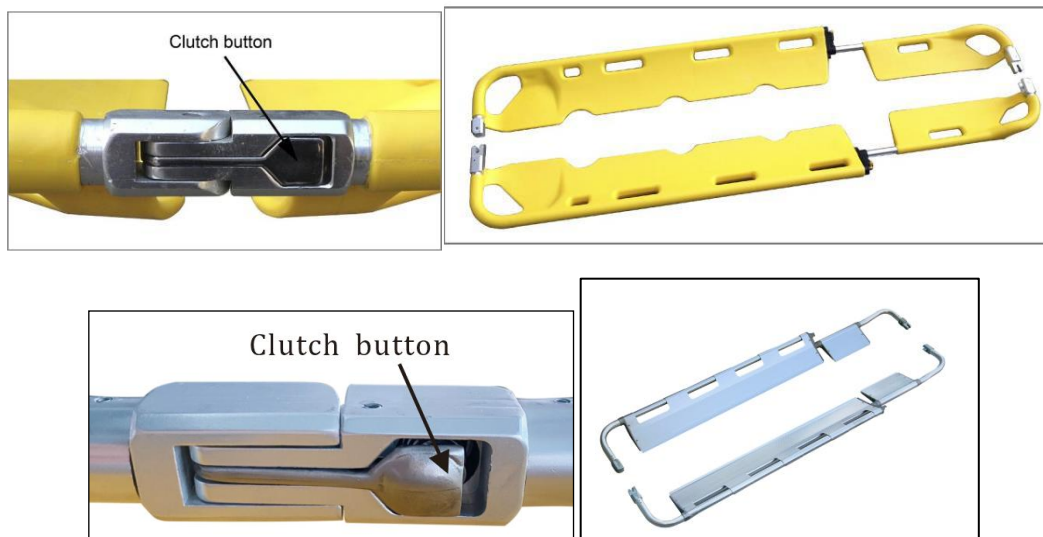Make sure the clutch button is completely locked when transporting patients.

Patient positioning:

Apply the stretcher following the minimal-handling principle ("skin to scoop"). Before moving the patient, an operator manually stabilizes the head in a neutral position and maintains this stabilization throughout the manoeuvre.

05 OPERATION

REV. 2026/05

1. Two operators position themselves on one side of the patient and, with minimal lateralization (no more than 15–20°), gently lift the patient's side.
2. A third operator carefully slides one half of the stretcher under the patient, against the back, without displacing the patient longitudinally.
3. Return the patient to the supine position and repeat the operation on the opposite side with the second half.
4. Bring the two halves together, aligning the upper and lower anchoring points. Engage the lower latch (foot end) first and then the upper latch (head end).
5. Before proceeding, visually check that the patient's skin, clothing, straps and life-support equipment have not been trapped between the two halves.
6. Verify that both latches are correctly engaged by visual inspection and by applying gentle traction at both ends.

Patient immobilization:

Apply the restraint straps supplied by the manufacturer (at least three) or, alternatively, a spider-strap system, distributing them across the chest, pelvis and thighs. Tighten the straps firmly, without compromising the patient's breathing and without compressing bony prominences or wounds.

The medical professional in charge shall decide, according to the patient's clinical condition, whether to apply a cervical collar, head immobilizer or ankle immobilizer as a supplement. If no head immobilizer is used, an operator must maintain manual cervical stabilization throughout the transfer.

Lifting and transfer:

Lifting is performed manually only. Never use cranes, slings, hoisting systems, rope rescue systems or in high-angle or vertical rescues.

The number of operators is determined by the weight of the patient, of the device itself and of the accessories. The experience gathered by Promeba advises performing the lift with four (4) operators whenever possible; in some cases, two (2) may be sufficient.

When only two operators are involved, they should be positioned centrally or slightly shifted towards the head end, where the greatest weight is concentrated. The placement must be symmetrical and allow safe, stable lifting of both the head and foot ends of the stretcher.

During lifting, keep the stretcher horizontal with respect to the ground. Abnormal tilts may cause harm to the patient, the device and the operators.


The maximum load per operator must never exceed the limit established by the local occupational health and safety regulations.

Once the main ambulance stretcher or vacuum mattress has been reached, transfer the patient and remove the scoop stretcher, as it is not suitable for sanitary transport in an ambulance.

Storage of the stretcher:

Once use is complete, the stretcher can be folded for storage in the corresponding compartment of the ambulance or in its usual location.

1. Unlock the adjustment buttons on both sides of the stretcher and hold them in the vertical position.
2. Extend the stretcher to its maximum length, ensuring that the extension joints are fully exposed.
3. Fold the head end over the foot end, considerably reducing the dimensions of the device.

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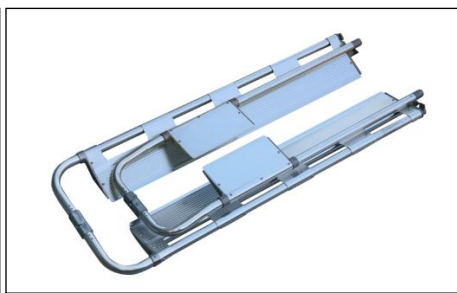06 GENERAL MAINTENANCE *REV. 2026/05*

06.1. *Cleaning*

Keeping the equipment clean is essential to ensure proper use and durability of the assembly. Thorough cleaning must be carried out periodically, especially in areas exposed to dirt such as the handles.

Do not use high-pressure cleaning systems or bleach to disinfect the product, as they may damage it. Instead, use an alcohol-based disinfectant, wash with water and allow it to dry naturally; do not use direct heat sources for drying.

Failure to perform cleaning operations may lead to the risk of cross-infection due to the presence of secretions and/or residues. During all inspection and disinfection operations, the operator must wear appropriate personal protective equipment such as gloves, eye protection, etc.

When the cleaning procedure includes a disinfectant product, carefully follow the instructions of the manufacturer of the product used regarding the method of application and contact time. Make sure that all necessary precautions have been taken to ensure that there is no risk of cross-infection or contamination of patients and operators.

06.2. *Lubricating*

In general, all moving parts must be greased. Our products leave the factory fully greased and lubricated, however, it is possible that over time and use of the product, the elements become degreased, either due to loss of lubrication or dirt.

Periodically clean and grease the affected areas according to the manufacturer's specifications. Check that there are no loose, missing or worn parts.

Periodically inspect all moving elements to ensure that the components are tight.

06.3. *Wear areas*

Periodic inspection of the different components of the device for signs of wear is a preventive measure that can reduce breakdowns. Pay particular attention to the following items as the main wear points: the restraint straps; the length-adjustment buttons; the upper and lower latches that join the two halves of the stretcher; the sliding extension profile; and the aluminium structure (checking for deformations, cracks or abrasions).

In the event that a worn component needs to be changed, you should contact our sales department for more information on ordering spare parts and their installation.

The repair or replacement of product components must necessarily be carried out by qualified personnel, that is, who have received training from the Manufacturer, Promeba S.L. and original spare parts should always be used.

Promeba S.L. declines any responsibility for any damage, direct or indirect, resulting from improper use of spare parts and/or any repair intervention carried out by unauthorized persons.

06 GENERAL MAINTENANCE *REV. 2026/05*

06.4. Maintenance program

Although the frequency of the controls is determined by factors such as legal requirements, the type of use, the frequency of use, the environmental conditions during use and storage, it is placed a recommendation on the maintenance that the equipment must have for correct operation throughout its useful life.

All maintenance and revision activities must be recorded and documented with the corresponding technical intervention reports, and the documentation must be kept for at least 10 years from the end of the product's useful life and must be made available to the competent Authorities and/or from the Manufacturer, when requested.

Maintenance review	Every use	When necessary	Every month	Every year
Disinfect	X			
Clean		X	X	X
Inspect	X	X	X	X

For both the disinfection and cleaning procedures, point 6.1 of this manual must be followed.

During the inspection procedure, which must be done after each use, the following tasks must be carried out:

- Check that all components are present.
- Check the integrity of the device, that is, that there are no breaks, cracks, holes or cuts
- Check the wear of the product following section 6.3 of this user guide.
- Check that the moving parts slide correctly.
- Check that the accessories provided are present, intact and that they function correctly.

06.5. Troubleshooting

Problem	Cause	Solution
The length-adjustment button does not lock the profile	Dirt or worn spring	Clean thoroughly. If the problem persists, immediately remove the device from service and contact the manufacturer.
The upper or lower latch does not engage when joining the two halves	Structural deformation or dirt in the latch housing	Clean thoroughly. If the problem persists, remove the device from service and contact the manufacturer.
The two halves do not separate when the release buttons are pressed	Locking mechanism obstructed or insufficient lubrication	Clean and lubricate the mechanism. If the problem persists, contact the manufacturer.
Accidental release of one of the latches during use	Improper closure or damaged latch	Immediately remove the device from service and contact the manufacturer.
Excessive bending or deflection of the stretcher under load	Damaged structural components or load exceeding the maximum capacity	Immediately remove the device from service and contact the manufacturer.
The extension profile does not slide or slides with difficulty	Dirt or lack of lubrication on the sliding surfaces	Clean and lubricate according to section 6.2 of this manual.
The stretcher cannot be folded for storage	The extension joints have not been fully exposed before folding	Extend the stretcher further until the joints are completely exposed, then fold.
Visible cracks, deformations or abrasions on the aluminium structure	Mechanical damage, impact or end of service life	Immediately remove the device from service and contact the manufacturer.
Frayed, cut or worn restraint straps	Normal wear, cuts or chemical degradation	Replace the straps with original spare parts from the manufacturer before the next use.

-

REV. 2026/05

Keep this document for at least 10 years after the end of the device's useful life.

51

08 MAINTENANCE REGISTER REV. 2026/05

Perform required maintenance as indicated by the manufacturer in this user manual.
Keep this document for at least 10 years after the end of the device's useful life.

DATE	SERVICE TYPE (Maintenance / Check / extension of life span)	OPERATIONS CARRIED OUT DURING THE MAINTENANCE	PERSON IN CHARGE OF THE SERVICE (Operator / Authorized / Centro / Manufacturer)

09 LEGAL NOTICES

REV. 2026/05

- This document may contain technical inaccuracies or typographical mistakes.
- Changes are periodically added to the information herein; these changes will be incorporated in new editions of the publication.
- Promeba, S.L. reserves the right to make any modification or improvement in the products described in this publication if it is appropriate.
- Promeba, S.L. may have patents or patent applications which address themes described in this document. The possession of this document does not entitle no license to these patents. The information contained in this document does not affect or change the specifications or product specifications or warranties of Promeba, S.L.
- No part of this document shall operate as express or implied license or compensation under the intellectual property rights of Promeba, S.L. or .third party.
- All the information contained in this document has been obtained in specific environments and it is presented as an illustration. The results obtained in other operating environments may vary.
- Promeba, S.L. can use or distribute the information that the customer provides in whatever way they see fit, without incurring any obligation with the client.

10 PRODUCT WARRANTY *REV. 2026/05*

Promeba, S.L. guarantees that its products have satisfactorily passed all established quality controls, both functional and material. The duration of the guarantee is 2 years from the date of purchase of the product.

This guarantee will be valid only when it is presented with the original invoice or proof of purchase (indicating the date of purchase, model and the name of the distributor) together with the defective product during the period covered by the guarantee. Promeba, S.L. reserves the right not to offer the free warranty service if the indicated documents are not presented or if the information they contain is incomplete or illegible.

1. This warranty will not apply if the model name or serial number of the product has been altered, erased, disappeared or becomes illegible.

2. This guarantee does not cover transportation costs or risks derived from the transportation of your product to and from Promeba, S.L.

3. This warranty does not cover any of the following:

- a) Periodic maintenance and repair or replacement of parts derived from normal wear and tear.
- b) Expendable material (components that are expected to need periodic replacement during the life of the product, such as belts, straps, etc.).
- c) Damages or defects derived from the improper use, operation or treatment of the product and not due to normal use of the product.
- d) Damages derived from:
 - i. Misuse, including:
 - Treatment that results in damages or physical, superficial or appearance changes of the product.
 - Installation, use or storage of the product in a way that does not

respect the instructions described by Promeba, S.L.

- Maintenance of the product in a way that does not respect the instructions of Promeba, S.L. for proper maintenance.

- Installation or use of the product in a way that does not respect the technical or safety regulations of the country where it is used or installed.

ii. Use of components not provided with the product or incorrect installation of accessory parts not previously tested.

iii. States or defects of the system in which the product is used or incorporated with the exception of other products Promeba, S.L. designed for use with the product.

iv. Use of the product with accessories, peripheral units and other products of a type, condition or standards not established by Promeba, S.L.

v. The manufacturer or distributor will be solely responsible for deciding whether to send the parts for repair, or to replace the product in its entirety. In no case will operators be sent to repair or replacement of the product.

Except in the cases mentioned above, Promeba, S.L. will make no warranties in relation to product, performance, accuracy, reliability, or fitness for logic purpose of the equipment or other type. If this exception is not legal or contemplated by current law, Promeba, S.L. will limit or exclude your warranties only to the extent permitted by applicable law.

The only obligation on the part of Promeba, S.L. in connection with this warranty is to repair or replace parts subject to the terms and conditions of this warranty.

Promeba, S.L. is not responsible for loss or damage to products, this warranty or others, including economic loss or non-assessable

10 PRODUCT WARRANTY *REV. 2026/05*

damage; the price paid for the product; loss of profits, income, information, usufruct or use of the product or associated products or indirect, accidental or critical loss or damage.

This clause refers to whether the loss or damage is due to deterioration or inoperability of the associated product due to defects or unavailability of Promeba, S.L., which has caused downtime, loss of user time or business interruption.

In cases where the law prohibits or limits these liability exclusions, Promeba, S.L. will exclude or limit his liability only to the extent permitted by applicable law.

For example, there are countries that prohibit the exclusion or limitation of damages caused by negligence, reckless negligence, willful misconduct, fraud and similar actions. The responsibility of Promeba, S.L. In this guarantee will not exceed, in any case, the price paid for the product, but if the current law allows only limitations of greater responsibilities, these will apply.

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*All rights reserved. Variations can be done whitout notice.
Promeba, S.L. is to be considered not responsible for damages caused by the lack or the wrongness of the information here mentioned.*

SCAN QR CODE

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Ctra C-16 Km 59.5 · 08650 Sallent (Barcelona) · SPAIN

T. 93 837 12 00

promeba@promeba.com · www.promeba.com